For Research Use Only

## CoraLite® Plus 488-conjugated Smac/DIABLO Polyclonal antibody



Catalog Number: CL488-10434 Featured Product

**Basic Information** 

Catalog Number: CL488-10434

Size: 1000 ug/ml Source: Rabbit Isotype:

diablo homolog (Drosophila)

Calculated MW: Immunogen Catalog Number: AG0696 27 kDa Observed MW:

20 kDa

BC011909

56616

Q9NR28

Full Name:

GeneID (NCBI):

**UNIPROT ID:** 

GenBank Accession Number:

**Purification Method:** 

Antigen affinity purification Recommended Dilutions: IF/ICC 1:50-1:500

Excitation/Emission maxima wavelengths:

493 nm / 522 nm

**Applications** 

**Tested Applications:** IF/ICC, FC (Intra) Species Specificity:

human, mouse, rat

Positive Controls: IF/ICC: HeLa cells,

## **Background Information**

Smac/DIABLO is a mitochondrial protein that promotes apoptosis by neutralizing members of the IAP family. SMAC/DIABLO is ubiquitously expressed in normal tissues, and the expression of SMAC messenger RNA (mRNA) in adult testis is the highest .Smac/DIABLO encodes a 239 amino acid protein (~27 kDa). The first 55 amino acids are a  $mit ochondrial\ targeting\ signal\ peptide\ that\ is\ cleaved\ to\ generate\ mature\ Smac-DIABLO\ (~21\ kDa).\ When\ released$ from mitochondria, Smac-DIABLO acts as a natural antagonist of inhibitor of apoptosis proteins (IAPs). Smac-DIABLO antagonistic activity is based on its N-terminal tetrapeptide (AVPI) that binds baculoviral IAP repeat (BIR) domains of IAPs, releasing their inhibitory effects on both initiator and effector caspases, thus promoting cell death. (PMID: 32325691, PMID: 32575872, PMID: 25650938)

Storage


Storage:

Store at -20°C. Avoid exposure to light. Stable for one year after shipment. Storage Buffer


PBS with 50% Glycerol, 0.05% Proclin300, 0.5% BSA, pH 7.3.

Aliquoting is unnecessary for -20°C storage

## Selected Validation Data

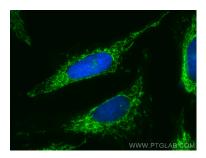

Immunofluorescent analysis of (4% PFA) fixed HeLa cells using CoraLite® Plus 488 DIABLO antibody (CL488-10434) at dilution of 1:100.

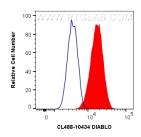

1X10^6 HeLa cells were intracellularly stained with 0.8 ug CoraLite® Plus 488 Anti-Human DIABLO (CL488-10434) (red), or 0.8 ug Isotype Control. Cells were fixed with 4% PFA and permeabilized with Flow Cytometry Perm Buffer (PF00011-C).