For Research Use Only

CoraLite® Plus 488-conjugated FABP5 Polyclonal antibody Catalog Number:CL488-12348 Featured Product


Basic Information	Catalog Number: CL488-12348	GenBank Accession Number: BC019385	Purification Method: Antigen affinity purification
	Concentration: 1000 ug/ml Source: Rabbit Isotype: IgG Immunogen Catalog Number: AG3005	GenelD (NCBI): 2171	Recommended Dilutions: IF/ICC: 1:50-1:500 FC (Intra): 0.40 ug per 10^6 cells in a 100 µl suspension Excitation/Emission maxima wavelengths: 493 nm / 522 nm
		UNIPROT ID: Q01469 Full Name: fatty acid binding protein 5 (psori associated) Calculated MW:	
		Applications	
IF/ICC, FC (Intra) Species Specificity: human, mouse, rat			
Background Information	FABP5, also named as PA-FABP and E-FABP, belongs to the calycin superfamily and Fatty-acid binding protein (FABP) family. It is high specificity for fatty acids. FABP5 is highest affinity for C18 chain length. It may be involved in keratinocyte differentiation. FABP5 is a fatty acid-binding protein and is expressed in epidermis and endothelial cells of the microvasculature of different organs. FABP5 has also been identified as a tumor-associated antigen, which is highly expressed in various cancers. FABP5 was detected in the sera of HNSCC patients with early stage cancer. Antibodies specific for FABP5 were significantly increased in a substantial amount in patients, suggesting that FABP5 may be a potential diagnostic biomarker for HNSCC. FABP5 may serve as a biomarker for HNSCC. (PMID:19602232)		
Storage	Storage: Store at -20°C. Avoid exposure to Storage Buffer: PBS with 50% glycerol, 0.05% Pro Aliquoting is unnecessary for -20°		nt.

For technical support and original validation data for this product please contact: T: 4006900926 E: Proteintech-CN@ptglab.com W: ptgcn.com This product is exclusively available under Proteintech Group brand and is not available to purchase from any other manufacturer.


proteintech

Antibodies | ELISA kits | Proteins www.ptglab.com

Selected Validation Data

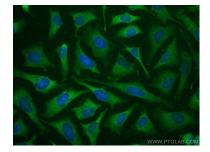

Immunofluorescent analysis of (-20°C Ethanol) fixed HeLa cells using Coralite® Plus 488 FABP5 antibody (CL488-12348) at dilution of 1:200.

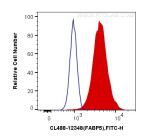

1X10^6 HeLa cells were intracellularly stained with 0.4 ug CoraLite® Plus 488 Anti-Human FABP5 (CL488-12348) (red), or 0.4 ug Control Antibody. Cells were fixed with 4% PFA and permeabilized with Flow Cytometry Perm Buffer (PF00011-C).