For Research Use Only

CoraLite® Plus 488-conjugated HDAC2 Polyclonal antibody



Catalog Number: CL488-12922

Featured Product

Basic Information

Catalog Number: CL488-12922

1000 µg/ml Source: Rabbit Isotype:

Immunogen Catalog Number:

AG3607

histone deacetylase 2 Calculated MW:

458 aa, 52 kDa; 488 aa,55 kDa

GenBank Accession Number:

BC031055

3066

Q92769

GeneID (NCBI):

UNIPROT ID:

Full Name:

Observed MW: 55-60 kDa

Purification Method:

Antigen affinity purification Recommended Dilutions: IF/ICC 1:50-1:500

Excitation/Emission maxima

wavelengths: 493 nm / 522 nm

Applications

Tested Applications: IF/ICC, FC (Intra) Species Specificity:

human, mouse, rat

Positive Controls:

IF/ICC: HepG2 cells,

Background Information

Histone deacetylases(HDAC) are a class of enzymes that remove the acetyl groups from the lysine residues leading to the formation of a condensed and transcriptionally silenced chromatin. Histone deacetylases act via the formation of large multiprotein complexes, and are responsible for the deacetylation of lysine residues at the N-terminal regions of core histones (H2A, H2B, H3 and H4). At least 4 classes of HDAC were identified. As a class I HDAC, HDAC2 was primarily found in the nucleus. HDAC2 forms transcriptional repressor complexes by associating with many different proteins, including YY1, a mammalian zinc-finger transcription factor. Thus, it plays an important role in transcriptional regulation, cell cycle progression and developmental events. This antibody is a rabbit polyclonal antibody raised against residues near the C terminus of human HDAC2.


Storage

Store at -20°C. Avoid exposure to light. Stable for one year after shipment. Storage Buffer


PBS with 50% Glycerol, 0.05% Proclin300, 0.5% BSA, pH 7.3.

Aliquoting is unnecessary for -20°C storage

Selected Validation Data

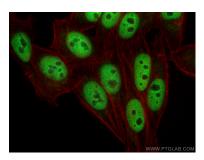

Immunofluorescent analysis of (4% PFA) fixed HepG2 cells using Coralite® Plus 488 HDAC2 antibody (CL488-12922) at dilution of 1:200, CL594-Phalloidin (red).

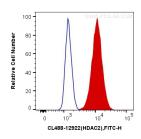

1X10^6 HEK-293T cells were intracellularly stained with 0.2 ug CoraLite® Plus 488 Anti-Human HDAC2 (CL488-12922) (red), or 0.2 ug Control Antibody. Cells were fixed and permeabilized with Transcription Factor Staining Buffer Kit (PF00011).