For Research Use Only

## CoraLite® Plus 488-conjugated iNOS Polyclonal antibody

Proteintech®
Antibodies | ELISA kits | Proteins
www.ptglab.com

Catalog Number: CL488-18985

Featured Product

**3 Publications** 

**Basic Information** 

Catalog Number:

CL488-18985

Size:

1000 µg/ml Source:

Rabbit Isotype:

løG

GenBank Accession Number:

NM\_000625

GeneID (NCBI):

4843

UNIPROT ID: P35228

Full Name:

nitric oxide synthase 2, inducible

Calculated MW: 131 kDa Observed MW:

110-130 kDa, 65-70 kDa

Purification Method:

Antigen affinity purification

Recommended Dilutions:

Excitation/Emission maxima wavelengths:

493 nm / 522 nm

IF/ICC 1:50-1:500

Applications

Tested Applications: IF/ICC, FC (Intra)

Cited Applications:

IF

Species Specificity: human, mouse, rat Cited Species:

human

Positive Controls:

IF/ICC: HepG2 cells,

## **Background Information**

NOS2, also named as iNOS and NOS2A, produces nitric oxide (NO) which is a messenger molecule with diverse functions throughout the body. NO is a reactive free radical which acts as a biologic mediator in several processes, including neurotransmission, antimicrobial and antitumoral activities. NOS2 is a nitric oxide synthase which is expressed in liver and is inducible by a combination of lipopolysaccharide and certain cytokines. iNOS has a very short half-life due to rapid degradation by calpain. iNOS monomer is a direct substrate of calpain I and can be cleaved by calpain I at the canonical CaM-binding site(503-532aa) of iNOS, and then a ~70kD band can be detected by western (PMID:11786228). This antibody is specific to NOS2.

## **Notable Publications**

Author	Pubmed ID	Journal	Application
Jiming Zhou	39185902	J Biochem Mol Toxicol	
Zeyan Liang	38488560	Neural Regen Res	IF
	38967729		

Storage

Storage


Store at -20°C. Avoid exposure to light. Stable for one year after shipment.

Storage Buffer:


PBS with 50% Glycerol, 0.05% Proclin300, 0.5% BSA, pH 7.3.

Aliquoting is unnecessary for -20°C storage

## Selected Validation Data

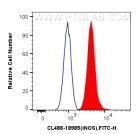

1X10^6 HepG2 cells were intracellularly stained with 0.4 ug CoraLite® Plus 488 Anti-Human iNOS (CL488-18985) (red), or 0.4 ug CoraLite® Plus 488-conjugated Rabbit IgG control Rabbit PolyAb (CL488-30000, Clone:) (blue). Cells were fixed with 4% PFA and permeabilized with Flow Cytometry Perm Buffer (PF00011-C).

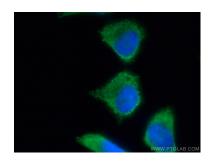

Immunofluorescent analysis of (-20°C Methanol) fixed HepG2 cells using CoraLite® Plus 488 iNOS antibody (CL488-18985) at dilution of 1:200.