

For Research Use Only



CoraLite® Plus 488-conjugated Phospho-Histone H2A.X (Ser139) Recombinant monoclonal antibody

Catalog Number:CL488-83307-2

Basic Information

Catalog Number:

CL488-83307-2

Source:

Rabbit

Isotype:

IgG

GenBank Accession Number:

BC013416

GeneID (NCBI):

3014

UNIPROT ID:

P16104

Full Name:

H2A histone family, member X

Calculated MW:

15 kDa

Observed MW:

15 kDa

Purification Method:

Protein A purification

CloneNo.:

5N19

Recommended Dilutions:

FC (Intra): 0.13 ug per 10⁶ cells in a
100 µl suspension

Excitation/Emission maxima
wavelengths:

493 nm / 522 nm

Applications

Tested Applications:

FC (Intra)

Species Specificity:

human

Positive Controls:

FC (Intra) : Staurosporine treated Jurkat cells,

Background Information

The histone variant H2AX is a major component of the DNA damage response (DDR), especially functioning in amplifying DNA damage signals. In response to DNA double-strand breaks (DSBs), H2AX is instantaneously phosphorylated at Ser139 (a form called γH2AX) by the kinases ATM and ATR. The phosphorylation of H2AX at Ser139, resulting in the formation of gamma-H2AX puncta in the nuclei, is an early event in the cellular response to DNA damage. Therefore, phospho-Histone H2A. X (Ser139) is also known as γ H2AX. The phosphorylation site of H2AX, Ser139, has also been described as Ser140 in other literature, and they recognize the same amino acid site. (PMID: 22908299, PMID: 30106130, PMID:22941631)

Storage

Storage:

Store at -20°C. Avoid exposure to light. Stable for one year after shipment.

Storage Buffer:

PBS with 50% glycerol, 0.05% Proclin300, 0.5% BSA, pH7.3

Aliquoting is unnecessary for -20°C storage

For technical support and original validation data for this product please contact:


T: 4006900926

E: Proteintech-CN@ptglab.com

W: ptgcn.com

This product is exclusively available under Proteintech Group brand and is not available to purchase from any other manufacturer.

Selected Validation Data

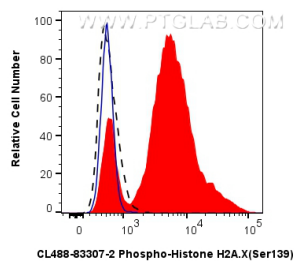

1x10⁶ Jurkat cells untreated (dashed lines) or treated with Staurosporine which intracellularly stained with 0.13 ug Coralite® Plus 488 Phospho-Histone H2a.X (Ser139) Recombinant Antibody (CL488-83307-2, Clone:5N19)(red), or 0.13 ug Coralite® Plus 488 Rabbit IgG Isotype Control RecAb (CL488-98136, Clone: 240953C9) (blue). Cells were fixed and permeabilized with Flow Cytometry Phosphorylated Protein Fixation/Permeabilization Kit (PF00026).