## For Research Use Only

## SUMO1-Trap Agarose, Kit for Immunoprecipitation



www.ptgcn.com

Catalog Number: s1tak

**Catalog Number: Basic Information** 

Alpaca **Applications:** IP, Co-IP Type: Nanobody Conjugate: Agarose beads; ~90 um (cross-linked 4% agarose beads) Class: Recombinant

The ChromoTek SUMO1-Trap Agarose, Kit for Immunoprecipitation consists of an anti-SUMO-Tag NANOBODY®/VHH, which is coupled to agarose beads. It also contains lysis, wash, and elution buffers that can be used for the immunoprecipitation of SUMO1 proteins from cell extracts of various organisms. **Description** 

**Host:** 

Binds specifically to human SUMO1 (Small Ubiquitin-like Modifier 1). The ChromoTek SUMO1-Trap Agarose will immunoprecipitate both monovalent SUMO1 and SUMO1 chains formed by the process of SUMOylation. In addition to human SUMO1, it will also bind SUMO1 homologs from other mammals. Specificity/Target

**Elution buffer** 2x SDS-sample buffer (Lammli), 200 mM glycine pH 2.5


Affinity (K<sub>D</sub>)

Storage

Storage: Shipped at ambient temperature. Upon receipt store at +4°C. Stable for one year. DO not freeze!

Storage Buffer: 20% ethanol

## **Selected Validation Data**

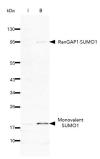

The SUMO 1-Trap Agarose (s1ta) was used to immunoprecipitate endogenous SUMO 1 and SUMO 1ylated proteins from HEK293T cells. Samples of the input lysate (I), and bound (B) fraction were analyzed by Western blot using 10329-1-AP (SUMO 1 Polyclonal Antibody). SUMO 1 and the commonly SUMO 1ylated protein RanGAP1 are highly enriched using the SUMO 1-Trap.